

Comment Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.

PCR protocol

Sheet ID PS_01818
Release date : 11.19.2008
Print date : 3/3/2014

RBRC 01386

Strain name Transgenic mouse, MuCre-A

Primer	primer_name	primer_seq	length		
1	Cre-ER F1	AACGAGTGATGAGGTTGCGAAGAA	24		Depositor's primer
2	Cre-ER R1	ATCAGCTACACCAGAGACGGAAA	23		Depositor's primer
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	Internal control (IL-2)	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	Internal control (IL-2)	JAX Genotyping Protocol


Reaction Components

Enzyme TaKaRa Ex Taq (Cat. No. RR001A)

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		15.3			
PCR buffer	10 X	2.5			
dNTP	2.5 mM	2.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
primer 4	10 pmol/ μ l	1.0			
DNA		1.0			
Taq		0.2			
total		25			

step	temp (°C)	time (sec)	note
1	94	180	
2	94	30*	
3	60	30*	
4	72	60*	*Repeat for 35 cycles
5	72	300	

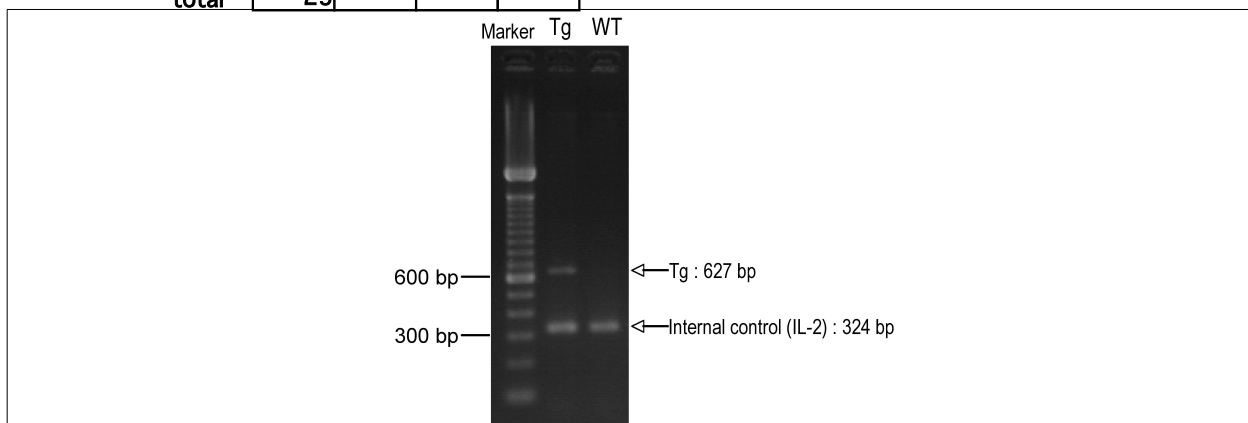

Marker

100 bp ladder (Invitrogen, Cat. 15628-050)

Gel

2% agarose in 1 X TAE (100 V, 35 min)

PCR products

primer_set	product_size	product
primer 1	627 bp	Tg
primer 3	324 bp	Internal control (IL-2)

Comment

PCR protocol

Sheet ID PS_00370
Release date : 11/22/2005
Print date : 3/3/2014

RBRC 01386

Strain name Transgenic mouse, MuCre-A

Primer	primer_name	primer_seq	length		
1	Cre-ER F1	AACGAGTGATGAGGTTGCGAAGAA	24		Depositor's primer
2	Cre-ER R1	ATCAGCTACACCAGAGACGGAAA	23		Depositor's primer
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	Internal control (IL-2)	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	Internal control (IL-2)	JAX Genotyping Protocol


Reaction Components

Enzyme Roche Taq (Cat. No. BM1418432)

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		16.3			
PCR buffer	10 X	2.5			
dNTP	2.5 mM	2.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	0.5			
primer 4	10 pmol/ μ l	0.5			
DNA		1.0			
Taq		0.2			
total		25			

step	temp (°C)	time (sec)	note
1	94	180	
2	94	30*	
3	60	30*	
4	72	60*	*Repeat for 30 cycles
5	72	300	

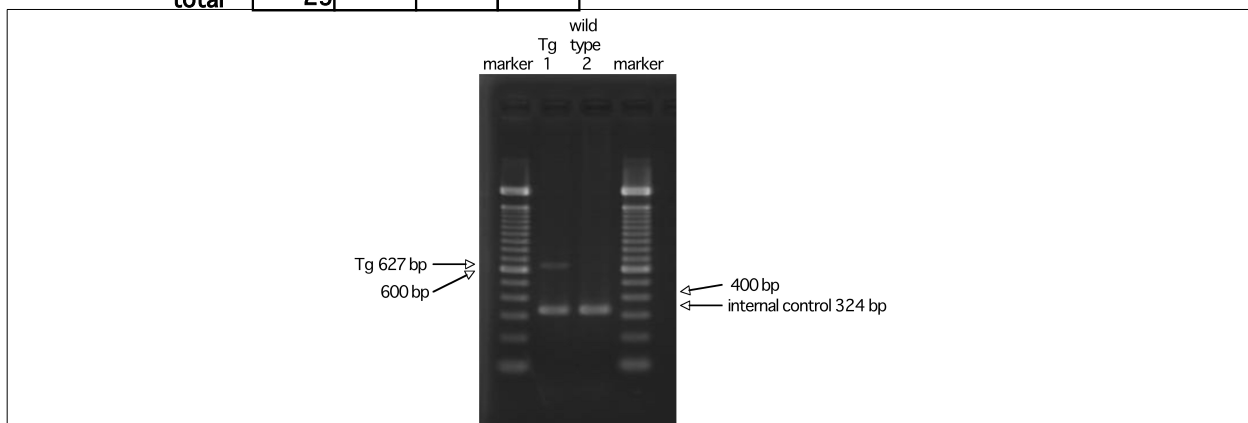

Marker

100 bp ladder (Invitrogen, Cat. 15628-050)

Gel

2% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set	product_size	product
primer 1	627 bp	Tg
primer 3	324 bp	Internal control (IL-2)

Comment