

PCR protocol

Sheet ID PS_05156
Release date : 03/21/2013
Print date : 3/21/2013

RBRC 02217

Strain name STOCK Tg(MoMuLV/RET-MEN2A)121Mat

Primer	primer_name	primer_seq	length		
1	HRET-1	CCTCCGCGAGAGCCGCAAAG	20		
2	94066	AATCCATGTGGAAGGGAGGGCTCGA	25		
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24		
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25		


Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		0.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
primer 4	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	1 U/ μ l	0.2			
total		20			

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

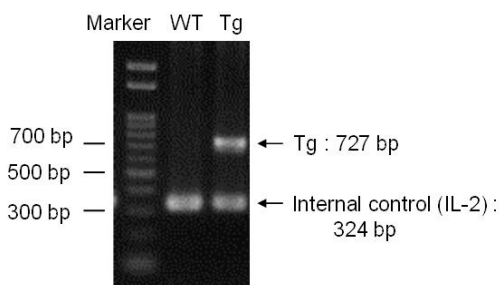

Marker

Gene Ladder 100 (Wako, Cat. 316-06951)

Gel

2% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set	product_size	product
primer 1	727 bp	Tg
primer 3	324 bp	Internal control (IL-2)

Comment