

PCR protocol

RBRC 03001

Strain name B6.Cg-Tg(CAG-INHBE)160hm

Primer	primer_name	primer_seq	length		
1	Sq-54	TGACCAAGCAGCAGCAAATCCTGGA	25		Depositor's primer
2	Sq-2	CTTTGAGGAGGCTGAAGACG	20		Depositor's primer
3	olMR0042	CTAGGCCACAGAATTGAAAGATCT	24		JAX Genotyping Protocol
4	olMR0043	GTAGGTGGAAATTCTAGCATCATCC	25		JAX Genotyping Protocol


Reaction Components

Enzyme TaKaRa Ex Taq (Cat. No. RR001A)

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		15.3			
PCR buffer	10 X	2.5			
dNTP	2.5 mM	2.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
primer 4	10 pmol/ μ l	1.0			
DNA		1.0			
Taq		0.2			
total		25			

step	temp (°C)	time (sec)	note
1	94	180	
2	94	30*	
2	60	30*	
4	72	60*	*Repeat for 30 cycles
5	72	300	

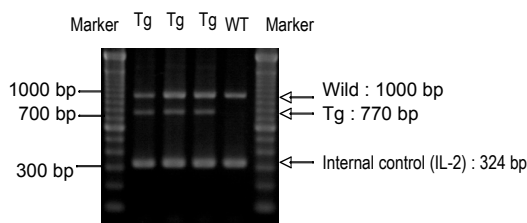

Marker

100 bp ladder (Invitrogen, Cat. 15628-050)

Gel

2% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set	product_size	product
primer 1	770 bp	Tg
primer 1	1000 bp	WT
primer 3	324 bp	Internal control (IL-2)

Comment