

PCR protocol

Sheet ID PS_05893
Release date : 7/11/2014
Print date : 7/11/2014

RBRC 03626

Strain name B6.B6D2-Tg(EF1A-BCL2)1Ymat

Primer	primer_name	primer_seq	length	
1	Ef1	GGTTTGCCGCCAGAACACAG	20	Depositor's primer
2	Bcl2	GGTTGACGCTCTCCACACAC	20	Depositor's primer
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	JAX Genotyping Protocol


Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		0.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
primer 4	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	1 U/ μ l	0.2			
total		20			

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	120*	*Repeat for 30 cycles

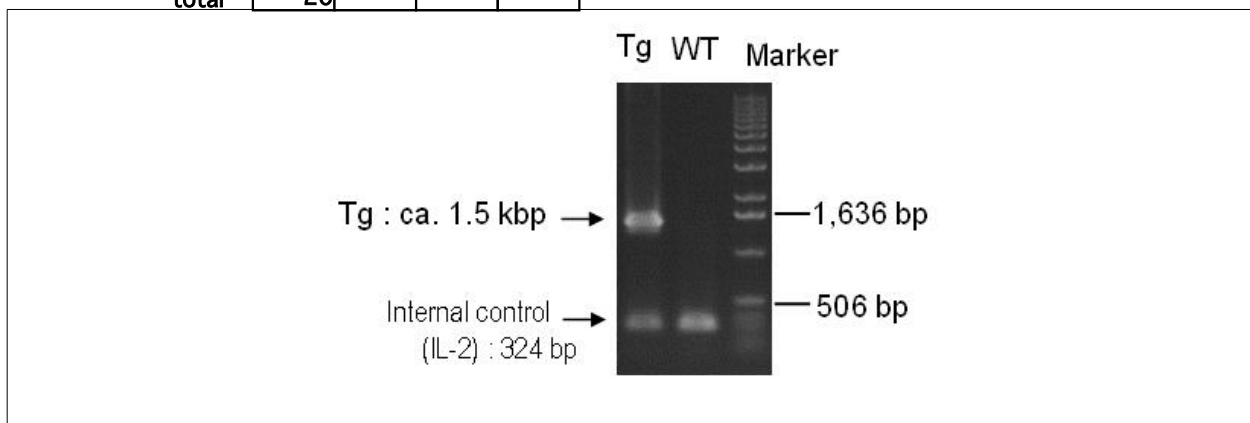

Marker

1 kb DNA ladder (Invitrogen, Cat. 15615-024)

Gel

1% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set	product_size	product
primer 1	primer 2	ca. 1.5 kbp
primer 3	primer 4	324 bp

Comment Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.

PCR protocol

Sheet ID PS_05593
Release date : 5/22/2014
Print date : 7/11/2014

RBRC 03626

Strain name B6.B6D2-Tg(EF1A-BCL2)1Ymat

Primer	primer_name	primer_seq	length	
1	Ef1	GGTTTGCCGCCAGAACACAG	20	Depositor's primer
2	Bcl2	GGTTGACGCTCTCCACACAC	20	Depositor's primer
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	JAX Genotyping Protocol


Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		0.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
primer 4	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	1 U/ μ l	0.2			
total		20			

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

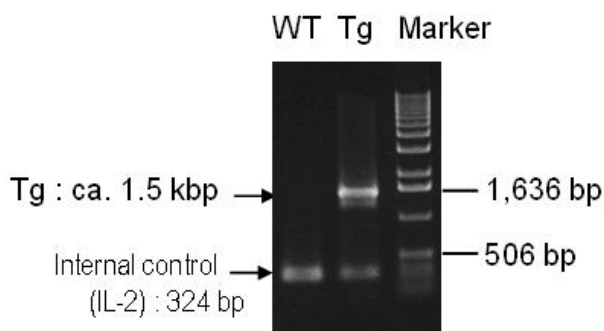

Marker

1 kb DNA ladder (Invitrogen, Cat. 15615-024)

Gel

1% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set	product_size	product
primer 1	primer 2	ca. 1.5 kbp
primer 3	primer 4	324 bp

Comment Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.