

RBRC 03830

Strain name STOCK Tg(Pitx2 17k-lacZ)26Hmd

Primer	primer_name	primer_seq	length	Reference
1	Pitx2c-S	GTTAGGCCAACAGGGAAGC	19	Depositor's primer
2	3' primer	GATGGGCGCATCGTAACCGTGC	22	Depositor's primer
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	JAX Genotyping Protocol

Reaction Components


Enzyme	TOYOBO KOD FX (Code.No : KFX-101)
--------	-----------------------------------

Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	90*	*Repeat for 30 cycles

Reagent	Conc	Reaction 1	2	3	4
H2O		0.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
primer 4	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	1 U/ μ l	0.2			
total		20			

Marker	2-Log DNA Ladder (NEB, Cat. N3200)
Gel	1% agarose in 1 X TAE (100 V, 30 min)

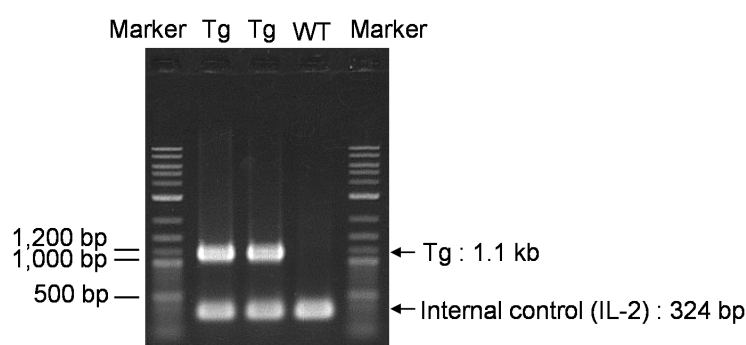

primer_set	product_size	product (allele type)
primer 1	primer 2	1.1 kb
primer 3	primer 4	324 bp

Non-Tg : 324 bp
Tg : 324 bp, 1.1 kb

Comment [Standard PCR](#)

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.

RBRC 03830

Strain name STOCK Tg(Pitx2 17k-lacZ)26Hmd

Primer	primer_name	primer_seq	length	Reference
1	ZS-6	GCACGGTTACGATGCGCCCCTACACCAACGT	33	Depositor's primer
2	ZA-6	ACGGCAAACGACTGTCCTGGCCGTAACCGACC	32	Depositor's primer
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	JAX Genotyping Protocol

Reaction Components


Enzyme	TaKaRa Ex Taq (Cat. No. RR001A)
--------	---------------------------------

Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	180	
2	94	30*	
3	60	30*	
4	72	60*	*Repeat for 30 cycles
5	72	300	

Reagent	Conc	Reaction 1	2	3	4
H2O		11.3			
PCR buffer	10 X	2.0			
dNTP	2.5 mM	1.6			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
primer 4	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	5U/ μ l	0.1			
total		20			

Marker	Gene Ladder 100 (Wako, Cat. 316-06951)
Gel	2% agarose in 1 X TAE (100 V, 30 min)

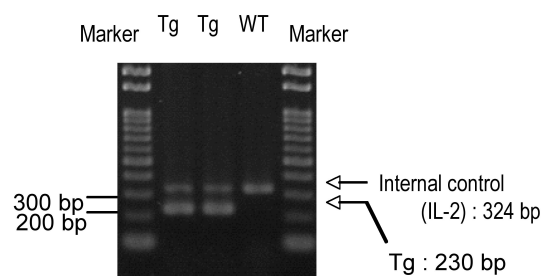

primer_set	product_size	product (allele type)
primer 1	primer 2	230 bp
primer 3	primer 4	324 bp

Non-Tg : 324 bp
Tg : 324 bp, 230 bp

Comment