

PCR protocol

RBRC 04017

Strain name C57BL/6J-Tg(CAG-floxed INHBE)80hm

Primer	primer_name	primer_seq	length	
1	Sq-54	TGACCAAGCAGCAGCAAATCCTGGA	25	Depositor's primer
2	Sq-2	CTTTGAGGAGGCTGAAGACG	20	Depositor's primer


Reaction Components

Enzyme	TaKaRa Ex Taq (Cat. No. RR001A)
--------	---------------------------------

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		13.3			
PCR buffer	10 X	2.0			
dNTP	2.5 mM	1.6			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	5U/ μ l	0.1			
total		20			

step	temp (°C)	time (sec)	note
1	94	180	
2	94	30*	
3	60	30*	
4	72	60*	*Repeat for 30 cycles
5	72	300	

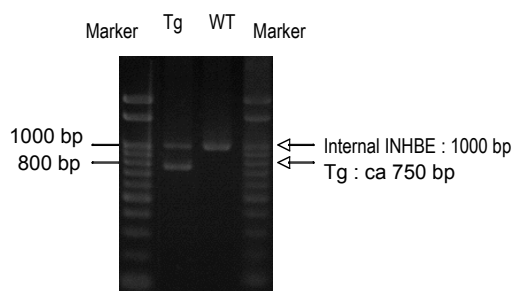

Marker

Gene Ladder 100 (Wako, Cat. 316-06951)

Gel

2% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set		product_size	product
primer 1	primer 2	ca 750 bp	Tg
primer 1	primer 2	1 000 bp	Internal INHBE

Comment