

RBRC 04110

Strain name B6;129-Cyp26a1<tm1.2Hmd>

Primer	primer_name	primer_seq	length	Reference
1	5C(CTV-5C)	TCCATGAGCCGAGTTCTGTAG	21	Genes Dev., 15, 213-225 (2001).
2	3A(CTV-3A)	TCCGTTTCAGTCAGTCTCATAC	21	Genes Dev., 15, 213-225 (2001).
3	5B(CTV-5B)	GGTAACGTGGGCAGTAACCTG	21	Genes Dev., 15, 213-225 (2001).

Reaction Components


Enzyme	TOYOBO KOD FX (Code.No : KFX-101)
--------	-----------------------------------

Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

Reagent	Conc	Reaction 1	2	3	4
H2O		1.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/μl	1.0			
primer 2	10 pmol/μl	1.0			
primer 3	10 pmol/μl	1.0			
DNA		1.0			
Taq	1 U/μl	0.2			
total		20			

Marker	Gene Ladder 100 (Wako, Cat. 316-06951)
Gel	2% agarose in 1 X TAE (100 V, 30 min)

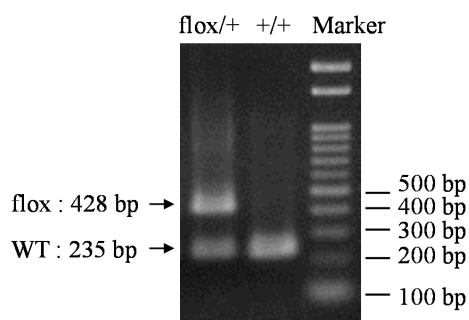

primer_set	product_size	product (allele type)
primer 1	primer 2	235 bp
primer 1	primer 2	428 bp
primer 2	primer 3	371 bp

Comment [Standard PCR](#)

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.