

PCR protocol

Sheet ID PS_06931
Release date : 5/11/2017
Print date : 5/11/2017

RBRC 04377

Strain name C57BL/6N-Zkscan3<Gt(LTOCG)050sb>

Primer	primer_name	primer_seq	length	
1		AGAATGTGTCGAGCAGAGTGGTG	23	RIKEN BRC
2	ARC07-02-08	TCTCTTCTGCTGCTCAGCCTTTC	23	Depositor's primer
3	4135	GATCTCTAGTTACCAGAGTCAC	22	Depositor's primer

Reaction Components


Enzyme	TOYOBO KOD FX (Code.No : KFX-101)
--------	-----------------------------------

Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		1.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/μl	1.0			
primer 2	10 pmol/μl	1.0			
primer 3	10 pmol/μl	1.0			
DNA		1.0			
Taq	1 U/μl	0.2			
total		20			

Standard PCR

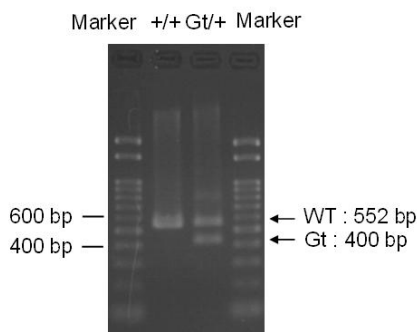

Marker

Gene Ladder 100 (Wako, Cat. 316-06951)

Gel

2% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set	product_size	product
primer 1	552 bp	WT
primer 3	400 bp	Gt

Comment

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.