

RBRC 04847

Strain name B6.B6129-Gt(ROSA)26Sor<tm1 (CAG-kikGR)Kqwa>

[illegible]

Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions


step	temp (°C)	time (sec)	note
------	--------------	---------------	------

[illegible]

Reagent	Conc	Reaction 1	2	3	4
H ₂ O		1.8			(μL)
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/ μl	1.0			
primer 2	10 pmol/ μl	1.0			
primer 3	10 pmol/ μl	1.0			
DNA		1.0			
Taq	1 U/ μl	0.2			
					Mari
					Gel
total		20			

Marker Gene Ladder 100 (Wako, Cat. 316-06951)

Gel	2% agarose in 1 X TAE (100 V, 30 min)
------------	---------------------------------------

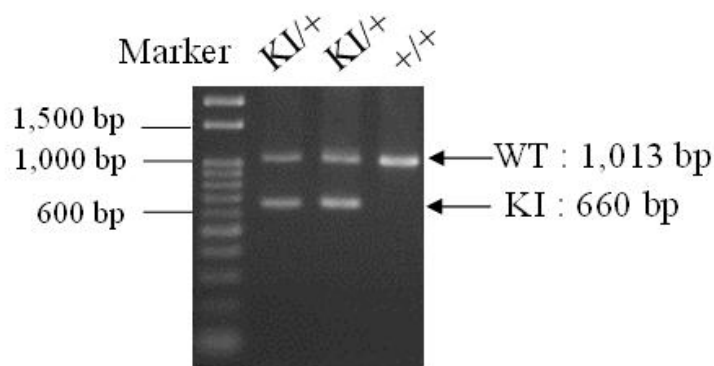

primer_set		product_size	product (allele type)
primer 1	primer 3	1,013 bp	WT
primer 1	primer 2	660 bp	KI

Comment Standard PCR

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.
Corrected primer 2 and 3 in the primer set on May 9, 2023.