

PCR protocol

Sheet ID PS_04128
Release date : 09/13/2013
Print date : 9/13/2013

RBRC 04900

Strain name B6;B6C3(129)-Tg(GATA6-mTomato)3Hmd

Primer	primer_name	primer_seq	length	
1		ACAGCGCTGGATGATTATGGGTCG	24	RIKEN BRC
2	Tomato-3'-2	TGAACTCCACCAGGTAGTGGC	21	Depositor's primer
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	JAX Genotyping Protocol


Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		0.8			
PCR buffer	2 X	10.0			
dNTP	2 mM	4.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
primer 4	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	1 U/ μ l	0.2			
total		20			

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

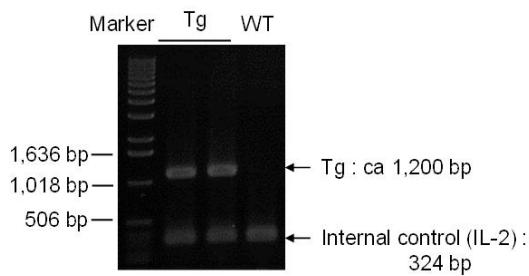

Marker

1 kb DNA ladder (Invitrogen, Cat. 15615-024)

Gel

1% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set	product_size	product
primer 1	primer 2	ca 1200 bp
primer 3	primer 4	324 bp

Comment