

PCR protocol

Sheet ID PS_04994
Release date : 9/12/2014
Print date : 9/12/2014

RBRC 05449

Strain name B6(129)-Tg(Plp1-tTS)1Kftnk

Primer	primer_name	primer_seq	length	
1	PLPU-604	TTTCCCATGGTCTCCCTTGAGCTT	24	Depositor's primer
2	mtTA24L	CGGAGTTGATCACCTTGGACTTGT	24	Depositor's primer
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	JAX Genotyping Protocol


Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		0.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
primer 4	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	1 U/ μ l	0.2			
total		20			

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

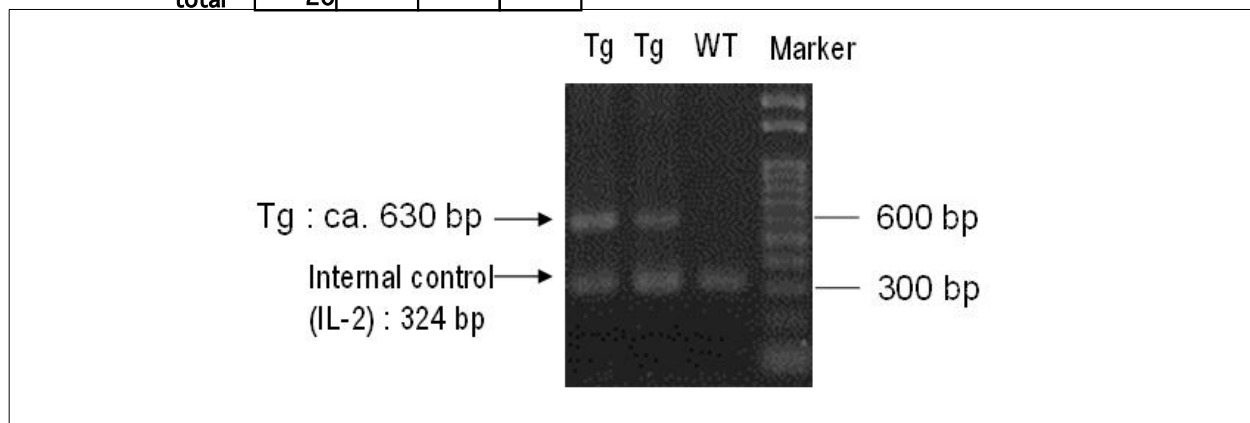

Marker

Gene Ladder 100 (Wako, Cat. 316-06951)

Gel

2% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set	product_size	product
primer 1	ca. 630 bp	Tg
primer 3	324 bp	Internal control (IL-2)

Comment Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.