

PCR protocol

Sheet ID PS_04765
Release date : 03.21.2013
Print date : 7/31/2014

RBRC 05451

Strain name B6;129-Nlgn3<tm1(STOP-tetO)Rhn>

Primer	primer_name	primer_seq	length	
1	NlgnU-465	GGGCTGGATGTTGCAATTGGAGTT	24	Depositor's primer
2	NlgnL420	TCCGTGGGCACATACACATTTCAGA	24	Depositor's primer
3	PGKF	GCGCCTACCGGTGGATGTGGAATGTGTGCG	30	


Reaction Components

Enzyme TaKaRa Ex Taq (Cat. No. RR001A)

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		12.3			
PCR buffer	10 X	2.0			
dNTP	2.5 mM	1.6			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
Taq		0.1			
DNA	5 U/ μ l	1.0			
total		20			

step	temp (°C)	time (sec)	note
1	94	180	
2	94	30*	
3	60	30*	
4	72	60*	*Repeat for 30 cycles
5	72	300	

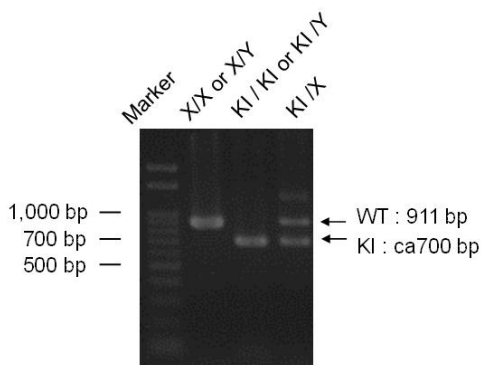

Marker

Gene Ladder 100 (Wako, Cat. 316-06951)

Gel

2% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set		product_size	product
primer 1	primer 2	911 bp	WT
primer 1	primer 3	ca 700 bp	KI

Comment