

PCR protocol

Sheet ID PS_06301
Release date : 8/14/2015
Print date : 8/14/2015

RBRC 06806

Strain name C57BL/6-Tg(hOX-EGFP,-cre)1Ahky

Primer	primer_name	primer_seq	length	
1	hOX-promoter_F	GCAGCGGCCATTCTTGG	18	Depositor's primer
2	Cre	CTAATCGCCATCTTCCAGCAGG	22	
3	mLC3ex3GT	TGAGCGAGCTCATCAAGATAATCAGGT	27	Mol. Biol. Cell., 15, 1101-1111(2004).
4	mLC3ex4AG	GTTAGCATTGAGCTGCAAGCGCCGTCT	27	Mol. Biol. Cell., 15, 1101-1111(2004).

Reaction Components


Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		0.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
primer 4	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	1 U/ μ l	0.2			
total		20			

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	180*	*Repeat for 30 cycles

Standard PCR

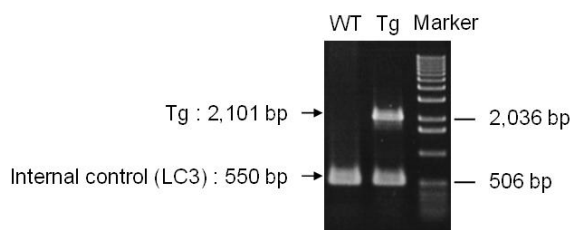

Marker

1 kb DNA ladder (Invitrogen, Cat. 15615-024)

Gel

1% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set	product_size	product
primer 1	primer 2	2,101 bp
primer 3	primer 4	550 bp

Comment Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.