



Sheet ID PS_07237
Release date : 7/29/2021
Print date : 7/29/2021


Strain name C57BL/6-Drd2<tm1.1Mmsh>

[illegible]

Enzyme	TOYOBO KOD FX (Code.No : KFX-101)
---------------	-----------------------------------

[illegible][illegible]

Gene Ladder 100 (Wako, Cat. 316-06951)
2% agarose in 1 X TAE (100 V, 30 min)

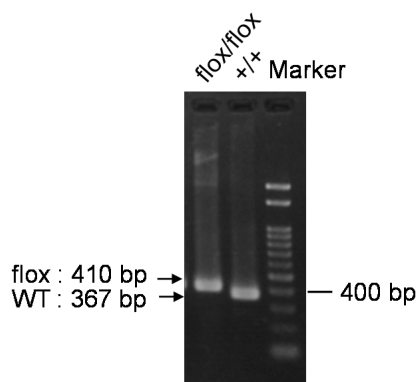

primer_set		product_size	product (allele type)
primer 1	primer 2	367 bp	WT
primer 1	primer 2	410 bp	flox

Comment Standard PCR

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.