

RBRC 10992

Strain name B6.Cg-Tg(Prnp-TARDBP*R208)#Koya

Primer	primer_name	primer_seq	length	Reference
1	promoter	GCCAGTCCGGAGCGGTGACTCA	22	RIKEN BRC
2	hTDP43-1012Rev	CTGCTATTACCAAATCCACCCTGATTC	27	Sci. Rep., 7(1):14972 (2017).
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	JAX Genotyping Protocol

Reaction Components

Enzyme	TOYOBO KOD FX (Code.No : KFX-101)
--------	-----------------------------------


Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	240*	*Repeat for 30 cycles

Reagent	Conc	Reaction 1	2	3	4
H2O		0.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
primer 4	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	1 U/ μ l	0.2			
total		20			

Marker	2-Log DNA Ladder (NEB, Cat. N3200)
Gel	1% agarose in 1 X TAE (100 V, 30 min)

Marker Tg WT Marker

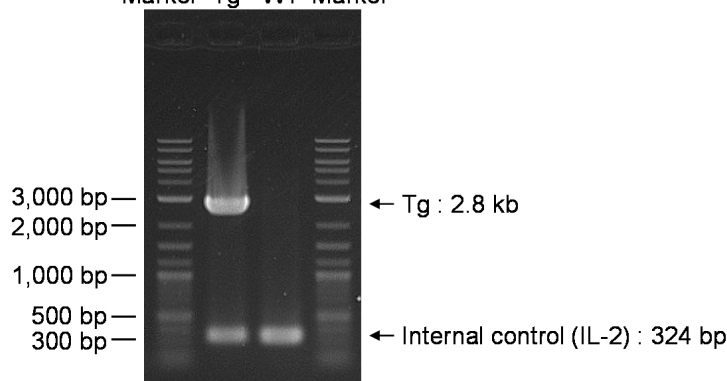

primer_set	product_size	product (allele type)
primer 1	primer 2	2.8 kb
primer 3	primer 4	324 bp

Non-Tg : 324 bp
Tg : 324 bp, 2.8 kb

Comment [Standard PCR](#)

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.