

RBRC 11302

Strain name C57BL/6N-Tg(CAG-Mpst)498Rmps

Primer	primer_name	primer_seq	length	Reference
1	pCAGGS_FW2	TAACCATGTTTCATGCCTTCT	20	EMBO Mol. Med., 11(12):e10695 (2019).
2	Mpst_Rv1	TTTACTGAGCCAGGGATGT	19	EMBO Mol. Med., 11(12):e10695 (2019).
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	JAX Genotyping Protocol

Reaction Components

Enzyme	TOYOBO KOD FX (Code.No : KFX-101)
--------	-----------------------------------


Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

Reagent	Conc	Reaction 1	2	3	4
H2O		0.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
primer 4	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	1 U/ μ l	0.2			
total		20			

Marker	Gene Ladder 100 (Wako, Cat. 316-06951)
Gel	2% agarose in 1 X TAE (100 V, 30 min)

Marker Tg WT Marker

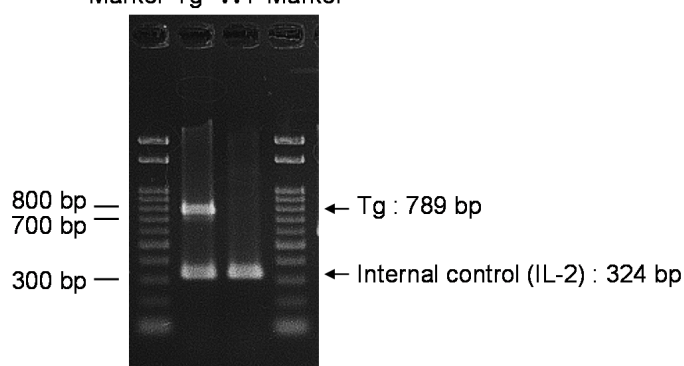

primer_set	product_size	product (allele type)
primer 1	primer 2	789 bp
primer 3	primer 4	324 bp

Non-Tg : 324 bp
Tg : 324 bp, 789 bp

Comment [Standard PCR](#)

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.